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
Nach EN 60601-1

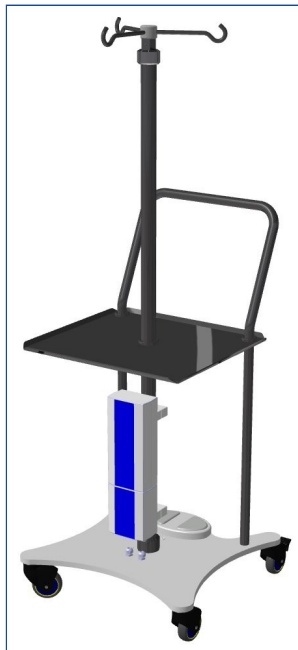Operating instructions

IV-Pole with electrical height adjustment

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1. General Information

Information relating to this manual

These operating instructions should enable you to use the IV-Pole with electrical height adjustment safely and properly. Please read these instructions carefully before using the device for the first time.

Pay particular attention to the instructions. If you give the device to a third party, please ensure that you also pass on the operating instructions at the same time.

Non-compliance with these instructions can lead to serious injury or may damage the device.

Correct use

The IV-Pole with electrical height adjustment is a lifting device for gravity infusions.

The device may not be operated in an explosive atmosphere or near strong magnetic fields.

Liability and warranty

All of the information and instructions provided are given to the best of our knowledge and are based on our experience and knowledge to date. We reserve the right to make technical changes essential for subsequent product development.

We refuse to accept any liability for injuries sustained as a result of non-compliance with these instructions.

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2. Symbols used in this operating instructions

Important instructions in this manual are highlighted with symbols, as shown below:

WARNING!

This symbol refers to safety instructions. Failure to comply with these could result in injury or even death.

WARNING! Electrical current!

Highlights hazardous situations caused by the electrical current. Failure to comply with this could result in injury or even death.

ATTENTION!

This refers to risks that may result in damage to the device.

NOTE!

Contains tips and other useful information listed in these operating instructions.

Disposal

When the life span of the device is over, especially when functional problems occur, dispose the device in accordance with environmental legislation applicable in your country. Electrical waste must not be disposed of with domestic waste. Contact your regional authorities or retailer for advice on recycling.

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3. Safety instructions:

WARNING!

It is absolutely necessary when setting up, handling and using the device that these instructions are followed precisely. If the system is inappropriately operated, it might cause danger for persons and objects!

- The devices are only to be used for the purposes named in the directions for use. Any other use may conceal unknown hazards and is to be refrained from.

WARNING! Electrical current!

- Any opening and exchange of electronic equipment may be carried out by appropriately trained personnel only.
- In case of malfunction, like e.g. development of smoke or gas, or unusual noises, all system components are to be turned off and all connections cables and mains cables are to be disconnected.
- The operating time as written in the Technical Data section may not be exceeded.
- Electric supply lines may not be exposed to danger of pinching, to bending or tensile stress.
- Modification of the control and the foot switch are not admissible.

WARNING!

- The device may not be operated in an explosive atmosphere!
- The allowed maximum weight per hook is 15 kg.
- The overall load must not exceed 30 kg.
- Lock all roller-brakes before lifting any load.
- The device must only be used in places where there is enough space above. Do not use the device under ceiling mounted hospital equipment.
- Make sure that there is enough space above the device, so that the pole can be lifted without colliding with any equipment mounted above. When operating the device always take care that no objects (e.g. furniture or body parts) get caught or jammed.
- Assurance against tipping (standing assurance) is guaranteed for up to 10% (5,7°) inclination for a load of 30 kg.

WARNING!

- Never lift the pole to maximum height if the inclination of the ground exceeds 5,7° and the maximum load is attached, because standing assurance can not be granted.
- When moving, make sure that the roller-brakes have been unlocked. Rolling across obstacles, soft flooring and slants must be done with particular care; it may also be necessary to lower the hanging load to the minimum infusion height.
- The electrically conductive rollers serve only to divert static charges.
- Provita GmbH is only responsible for the safety and proper functioning of the IV-Poles when the directions for assembly and use have been followed in all points and only when Provita accessories are used.

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4. Commissioning

4.1 Charging the battery

Connect the cord of the wall mounted battery charger to a suitable power outlet. The indicator “ON” is lit.

Connect the battery to charger.

The indicator “Charge” is lit.

When the battery is charged completely the indicator “Charge” turns off.

NOTE!

When the battery gets low an acoustic signal can be heard.

If so, disconnect the battery from the “Lifter” and recharge it with the battery charger.

4.2 Setup

- Position the device at the designated usage site and lock all roller-brakes.

WARNING!

Make sure that there is enough space above the device, so that the pole can be lifted without colliding with any equipment mounted above.

- IV-Poles not in use are always to be locked with the lock rollers.

4.3 Operation

4.3.1 Lifting a load

Use your foot to press the arrow up button (left) on the foot switch. Hold down the button until the Pole has reached the desired position.

NOTE!

The pole moves up until you release the switch, or until the maximum position is reached.

4.3.2 Sinking the load

Use your foot to press the arrow down button (right) on the foot switch. Hold down the button until the Pole has reached the desired position.

NOTE!

The pole moves down until you release the switch, or until the minimum position is reached.

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5. Box contents:

1x IV-Pole with electrical drive
1x Directions for assembly and use
1x rechargeable battery
1x wall mounted battery charger

6. Cleaning and Disinfection:

WARNING! Electrical current!

Always detach the plug from the power outlet before cleaning the device.

The infusion stand can be treated with mild cleaning agents, for the disinfection you should use commercially available disinfection agents for example: Amocid, Baktolan, Sagrotan.

7. Maintenance and Repair:

The infusion stand should be checked at least once a year for proper functionality of all important components:

To assure functionality of the locking and antistatic rollers clean them up.

Repairs have to be done by authorized personal or customer service with original equipment only.

8. Technical Details:

Input voltage: (Battery charger)	AC 230 V $\pm 6\%$ / 50 Hz
Operational voltage (electronics):	DC 24 V, 1x 2,9 Ah
Operating mode:	10% (1min on/ 9min off)
Protection class:	II
Type of protection:	IP 20
Temperature:	10 to 35°C
Moisture:	5 to 75% (not condensing)
Weight:	36,5 kg
Dimension (L,W,H1, H2):	~550,550,1500,2500 mm

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